


STAINLESS STEEL STERILIFE® SYRINGES



STERILIFE®

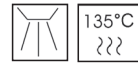

100% STAINLESS STEEL

1.8 cc Metric thread

1.8 cc American thread

2.2 cc American thread

| | 1.8 cc Metric thread | 1.8 cc American thread | 2.2 cc American thread |
|-----------------------------|----------------------|------------------------|------------------------|
| Réf. 2310 ¹ | • | | |
| Réf. 2310US ¹ | | • | |
| Réf. 2310-22US ¹ | | | • |
| Réf. 2320 ² | • | | |
| Réf. 2320US ² | | • | |
| Réf. 2320-22US ² | | | • |
| Réf. 2330 ³ | • | | |
| Réf. 2330US ³ | | • | |
| Réf. 2330-22US ³ | | | • |

¹ Lateral insertion of the cartridge without aspiration.


² Lateral insertion of the cartridge self aspiration.

³ Lateral insertion of the cartridge manual aspiration.

REF. 2310 : SYRINGE WITHOUT ASPIRATION

→ The piston rod exercises a constant pressure on the cartridge.

REF. 2320 : SELF-ASPIRATING SYRINGE

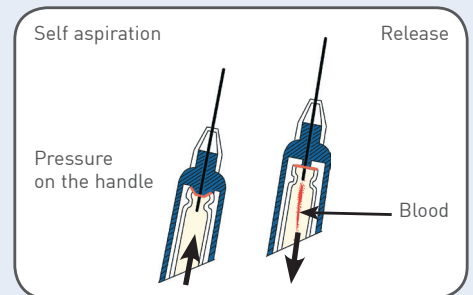

→ Thanks to a protrusion located at the bottom of the syringe tube, a pressure on the handle decreases the volume of anesthetic product in the cartridge.

→ Loosening the handle creates an air intake, and so an aspiration.

REF. 2330 : SYRINGE WITH MANUAL ASPIRATION

→ Manual aspiration thanks to a corkscrew system.